

CREATE A

# Signature LOOK

WITH THE WORLD'S  
**BROADEST** RANGE  
OF FILLERS

LIFT



FILL



VOLUMIZE



REFRESH



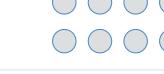
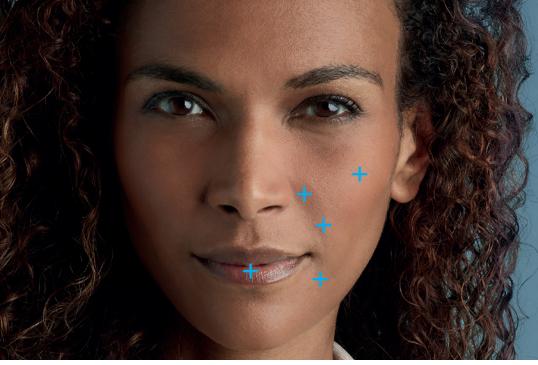
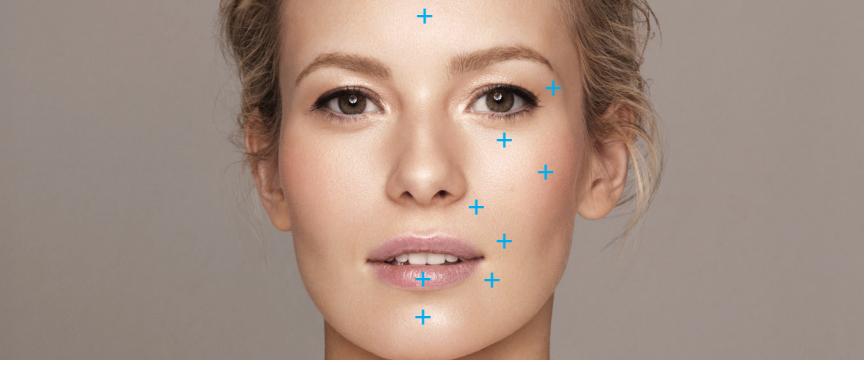
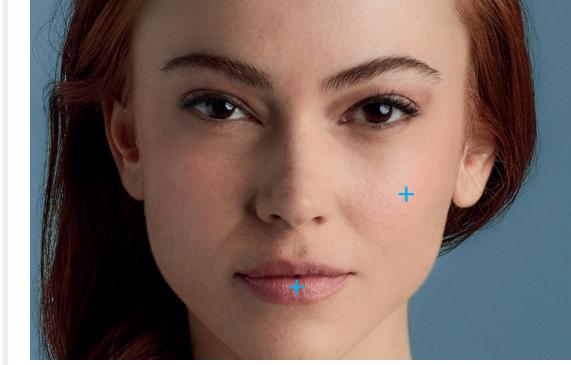
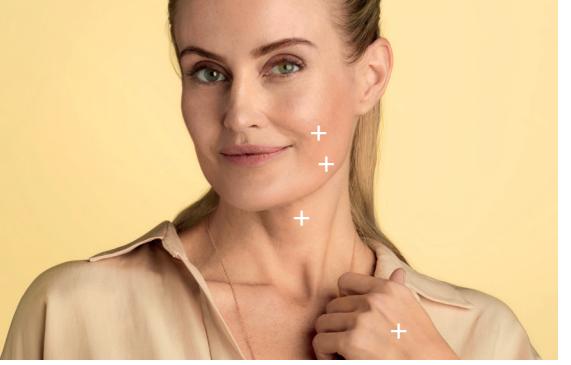
 GALDERMA

Restylane®

Restylane®  
SKINBOOSTERS™

# TWO PATENTED COMPLEMENTARY TECHNOLOGIES - NASHA™ AND OBT™ -

## CREATE A DIVERSE PORTFOLIO TO ACHIEVE A NATURAL LOOK

	LIFT		FILL			VOLUMIZE		REFRESH	
PRODUCT	Restylane Lyft Lidocaine <sup>1</sup> 	Restylane Defyne <sup>2</sup> 	Restylane Lidocaine <sup>3</sup> 	Restylane Silk <sup>4</sup> 	Restylane Refyne <sup>5</sup> 	Restylane Kysse <sup>6</sup> 	Restylane Volyme <sup>7</sup> 	Restylane Skinboosters Vital <sup>8</sup> 	Restylane Skinboosters Vital Light <sup>9</sup> 
SYRINGE VOLUME	1 mL	1 mL	0.5 & 1 mL	1 mL	1 mL	1 mL	1 mL	1 mL	1 mL
TECHNOLOGY	NASHA™	OBT™	NASHA™	NASHA™	OBT™	OBT™	OBT™	NASHA™	NASHA™
GEL FIRMNESS	Firm	Medium firm	Firm	Firm	Soft	Medium soft	Medium soft	Soft	Soft
GEL PARTICLE SIZE	 	 							
LIFTING CAPACITY	High	High	Moderate	Low	Moderate	Moderate	High; Volumizing	N/A	N/A
INTENDED USE & TISSUE TYPE									
	<ul style="list-style-type: none"> <li>For shaping facial contours, correction of folds &amp; lip enhancement</li> <li>For supportive tissue with severe lines and volume enhancement of cheeks</li> </ul>	<ul style="list-style-type: none"> <li>For correction of severe facial wrinkles &amp; folds such as nasolabial folds &amp; marionette lines</li> <li>For soft tissue with severe lines/deeper skin</li> </ul>	<ul style="list-style-type: none"> <li>For correction of wrinkles &amp; lip enhancement</li> <li>For supportive tissue with moderate lines</li> </ul>	<ul style="list-style-type: none"> <li>For restoring facial contours &amp; correction of thin superficial lines</li> </ul>	<ul style="list-style-type: none"> <li>For correction of moderate facial wrinkles &amp; folds such as nasolabial folds</li> <li>For soft tissue with moderate lines</li> </ul>	<ul style="list-style-type: none"> <li>For lip volume &amp; definition</li> </ul>	<ul style="list-style-type: none"> <li>For augmenting cheek volume</li> </ul>	<ul style="list-style-type: none"> <li>For improving skin smoothness, appearance &amp; elasticity in lower cheek/jawline &amp; dorsal hands</li> </ul>	<ul style="list-style-type: none"> <li>For improving skin elasticity in lower cheek/jawline &amp; upper neck</li> </ul>
INJECTION DEPTH	Deep dermis or surface of subcutis	Mid to deep dermis	Mid-dermis	Superficial dermis	Mid to deep dermis	Submucosal tissue	Supraperiosteic zone or subcutis	Dermal layer, preferably in the deeper part of the dermis	Dermal layer
NEEDLE RECOMMENDATION	29G <sup>1/2</sup> or 27G <sup>1/2</sup> (TWN)	27G <sup>1/2</sup> (UTWN)	29G <sup>1/2</sup> or 30G <sup>1/2</sup> (TWN)	29G <sup>1/2</sup> or 30G <sup>1/2</sup> (TWN)	30G <sup>1/2</sup> (UTWN)	30G <sup>1/2</sup> (UTWN)	27G <sup>1/2</sup> (UTWN)	29G (TWN)	29G (TWN)
MICROCANNULA RECOMMENDATION	23G-25G sterile blunt cannula	27G	27G-28G sterile blunt cannula	N/A	27G	25G-27G	27G	30G sterile blunt cannula	30G sterile blunt cannula

TWN=thin-walled needle; UTWN=ultra thin-walled needle.

NASHA™ technology creates a precise, firmer gel texture intended for a more pronounced fit.

OBT™ (Optimal Balance Technology™) creates a softer gel texture, allowing for smoother integration into the tissue for natural-looking results.

# TAILOR EVERY treatment



## IMPORTANT SAFETY INFORMATION FOR THE RESTYLANE RANGE OF FILLERS<sup>1-9</sup>

- Do not use on patients with known hypersensitivity to lidocaine or amide type local anaesthetics.
- Do not use on patients presenting with porphyria or bleeding disorders, or in patients who are taking thrombolytics or anticoagulants.
- Do not use if known history of keloids (for Restylane Skinboosters).
- Do not inject intravascularly or where an implant other than hyaluronic acid has been placed. Introduction of product into the vasculature may lead to embolization, occlusion of the vessels, ischemia, or infarction.
- Do not resterilize, overcorrect or mix with other products.
- Do not use in or near anatomic sites where there is active skin disease, inflammation, infection, tumours or related conditions.
- Rare but serious adverse events associated with the intravascular injection of soft tissue fillers in the face have been reported and include temporary or permanent vision impairment, blindness, cerebral ischemia or cerebral hemorrhage, leading to stroke, skin necrosis, and damage to underlying facial structures.
- Inject the product slowly and apply the least amount of pressure necessary. Immediately stop the injection if a patient exhibits any of the following symptoms: changes in vision, signs of a stroke, blanching of the skin, or unusual pain during or shortly after the procedure.
- Patients should receive prompt medical attention and possibly evaluation by an appropriate healthcare practitioner specialist should an intravascular injection occur.
- Common injection-related reactions such as erythema, swelling, pain, itching, discoloration, tenderness and temporary palpable lumpiness at the implant site may occur.
- In order to minimize the risks of potential complications, this product should only be used by healthcare practitioners who have appropriate training and experience, and who are knowledgeable about the anatomy at and around the site of injection.
- Healthcare practitioners are encouraged to discuss all potential risks of soft tissue injection with their patients prior to treatment and ensure that patients are aware of signs and symptoms of potential complications.

This is only a partial list of safety information. For complete product and safety information, please refer to the Instructions of Use.

## REFERENCES:

1. Restylane® Lyft Lidocaine Instructions for Use. **2.** Restylane® Defyne Instructions for Use. **3.** Restylane® Lidocaine Instructions for Use.
4. Restylane® Silk Instructions for Use. **5.** Restylane® Refyne Instructions for Use. **6.** Restylane® Kyssé Instructions for Use.
7. Restylane® Volume Instructions for Use. **8.** Restylane® Skinboosters™ Vital Lidocaine Instructions for Use.
9. Restylane® Skinboosters™ Vital Light Lidocaine Instructions for Use.



©2020 Galderma Canada Inc. All rights reserved. All trademarks are the property of their respective owners.  
RES-003-0120E

